




### 45mm x 2.4mm Full Ring - Mini

**Part No. 4.14501**

Expands capability of IMEX™ mini CESF system to include patients larger than appropriate for 35mm rings. Multiple hole design requires use of slotted/cannulated wire fixation bolts (part no. 4.30001 - shown below).

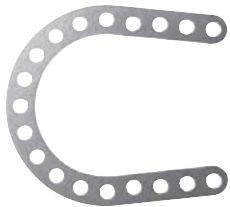

### 45mm x 2.4mm Stretch Ring - Mini

**Part No. 4.14504**


Expands capability of IMEX™ mini CESF system. Popular around hock and elbow to allow range of motion without a larger ring diameter being used. Multiple hole design requires use of slotted/cannulated wire fixation bolts (part no. 4.30001).



### Mini Angular Hinge Assembly

**Part No. 4.60001**

Can bolt to mini (CESF) rings to create strut and supports for hybrid frames. Also accepts mini hybrid rods in each end to create mini transarticular linear ESF device.

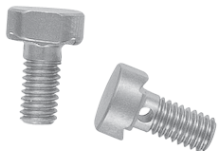

### Mini Slotted/Cannulated Wire Fixation Bolt

**Part No. 4.30001**

Bolt is slotted on each side and cannulated, providing three different points to capture and anchor fixation wires to rings. Designed to grip pins and wires from 0.8mm - 1.6mm with 0.8mm - 1.4mm being optimum. Required with new 45mm mini rings and optional with 35mm mini rings. Stainless steel with 4mm thread and 7mm wrench flats.



### Mini Slotted Wire Fixation Washer

**Part No. 4.30005**

Allows capture of fixation wires to rings (or rods) at locations otherwise occupied by supporting elements. Also utilized for drop wire technique. Can be pre-placed in several locations for optimal freedom. Designed to grip pins and wires from 0.8mm - 1.6mm with 0.8mm - 1.4mm being optimum. Multiple slots ensure washer is always in position for ease of use.



### Mini One-Hole Post

**Part No. 4.70005**

Mounts to mini rings for supporting additional wires, half-pins, and other elements. Stainless Steel.

**CALL TODAY to ORDER 800-828-4639**

1001 McKesson Drive • Longview, Texas 75604 • Phone 903-295-2196 • Fax 903-295-0730  
www.imexvet.com • imexdesk@imexvet.com